

## DDT Technical Mixture (500 ppm)

40 Years After The Ban: DDT-Related Compounds Accumulating in  
Southern California Bottlenose Dolphins

Susan A. Mackintosh, Nathan G. Dodder, Nellie J. Shaul, Lihini I. Aluwihare, Keith A. Maruya,  
Susan J. Chivers, Kerri Danil, David W. Weller, Eunha Hoh

Web Reference: <http://OrgMassSpec.github.io>

Prepared: 2015-03-16 12:28:50

SpecLibDDT version

OrgMassSpecR version 0.4-4

R version 3.1.0 (2014-04-10)

Name: chlorobenzaldehyde


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C<sub>7</sub>H<sub>5</sub>ClO

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. 2-chlorobenzaldehyde. All have the same mass spectrum, but exceed the number of possible isomers.

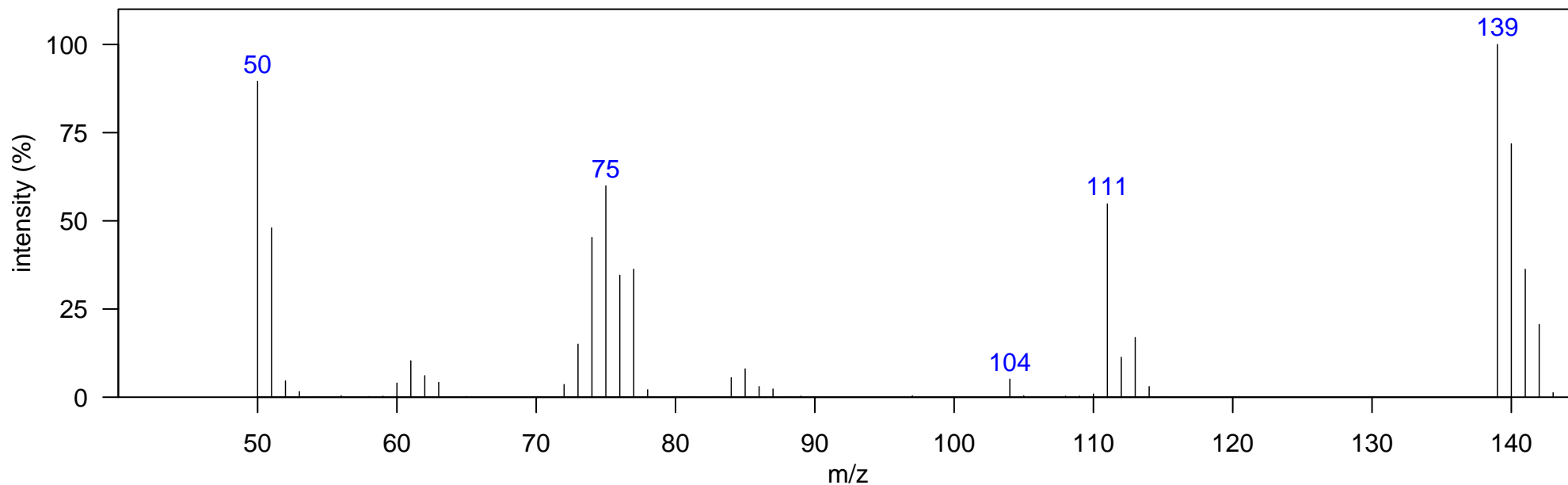

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	chlorobenzaldehyde	1	FALSE	FALSE	662.964 0.983
2	chlorobenzaldehyde	2	FALSE	FALSE	669.960 0.990
3	chlorobenzaldehyde	3	FALSE	FALSE	890.334 1.076
4	chlorobenzaldehyde	4	FALSE	FALSE	935.808 1.076
5	chlorobenzaldehyde	5	FALSE	FALSE	1086.220 0.766
6	chlorobenzaldehyde	6	FALSE	FALSE	1086.220 0.878
7	chlorobenzaldehyde	7	FALSE	FALSE	1110.710 0.898

m/z [Fragment]

75 [M-CHClO]<sup>+</sup>  
111 [M-CO]<sup>+</sup>  
139 M<sup>+</sup>

Name: trichloroethenylbenzene


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C<sub>8</sub>H<sub>5</sub>Cl<sub>3</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. trichloroethenylbenzene. All have the same mass spectrum, but exceed the number of possible isomers.

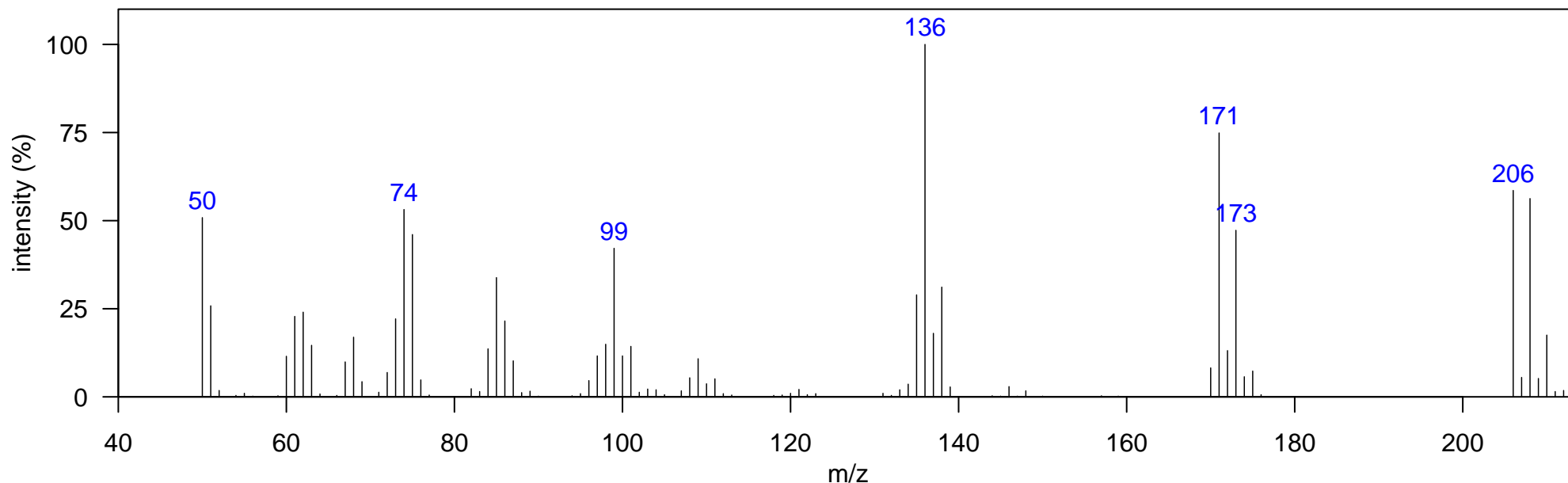

Isomer Information	Isomer InDolphin Monitored RT1D RT2D					
	1	trichloroethenylbenzene 1	FALSE	FALSE	872.844	0.990
	2	trichloroethenylbenzene 2	FALSE	FALSE	907.824	0.997
	3	trichloroethenylbenzene 3	FALSE	FALSE	918.318	1.003
	4	trichloroethenylbenzene 4	FALSE	FALSE	1054.740	0.845
	5	trichloroethenylbenzene 5	FALSE	FALSE	1093.220	0.858
	6	trichloroethenylbenzene 6	FALSE	FALSE	1436.020	0.713
	7	trichloroethenylbenzene 7	FALSE	FALSE	1600.430	0.700
	8	trichloroethenylbenzene 8	FALSE	FALSE	1600.430	1.696
	9	trichloroethenylbenzene 9	FALSE	FALSE	1610.920	1.676
	10	trichloroethenylbenzene 10	FALSE	FALSE	1614.420	0.713
	11	trichloroethenylbenzene 11	FALSE	FALSE	1635.410	0.700

m/z [Fragment]

136 [M-Cl<sub>2</sub>]<sup>+</sup>  
171 [M-Cl]<sup>+</sup>  
206 M<sup>+</sup>

Name: trichloro-(chlorophenyl) ethane


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C<sub>8</sub>H<sub>8</sub>Cl<sub>4</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. 1,2,2-Trichloro-1-(4-chlorophenyl)ethene

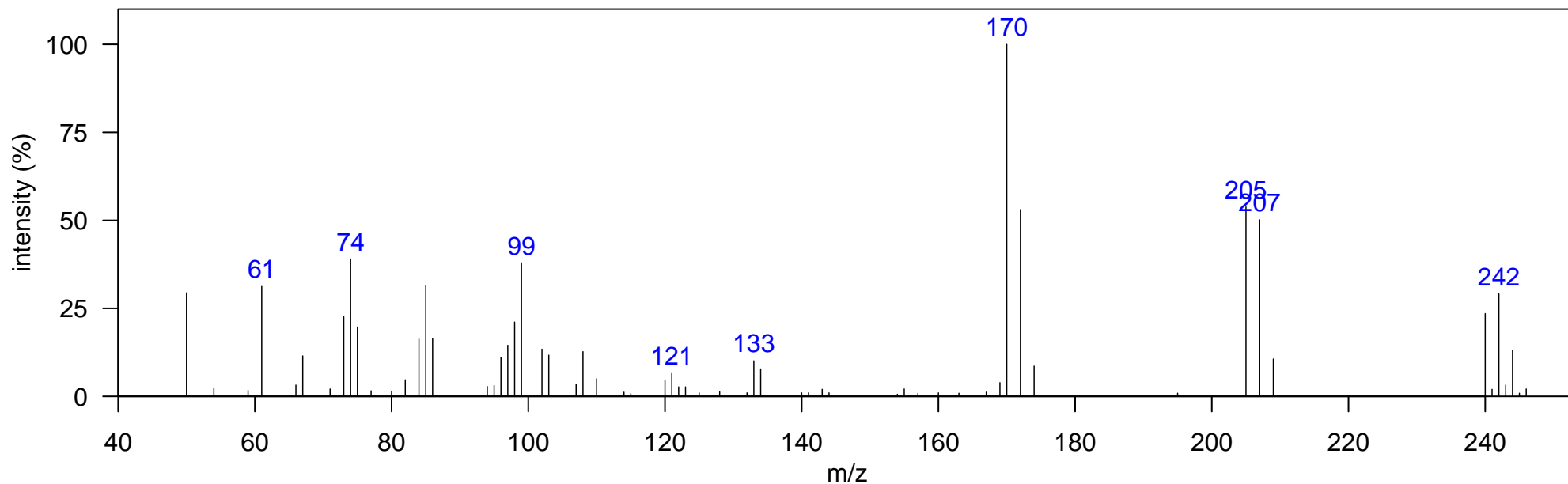

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	trichloro-(chlorophenyl) ethane 1	FALSE	FALSE	928.812	1.010
2	trichloro-(chlorophenyl) ethane 2	FALSE	FALSE	977.784	1.010
3	trichloro-(chlorophenyl) ethane 3	FALSE	FALSE	1652.900	1.676
4	trichloro-(chlorophenyl) ethane 4	FALSE	FALSE	1705.370	1.650

m/z [Fragment]

135 [M-Cl<sub>3</sub>]<sup>+</sup>  
170 [M-Cl<sub>2</sub>]<sup>+</sup>  
205 [M-Cl]<sup>+</sup>  
240 M<sup>+</sup>

Name: trichloro-(propenyl) benzene


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C<sub>9</sub>H<sub>7</sub>Cl<sub>3</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. Benzene, (3,3,3-trichloro-1-propenyl)-, (E)-

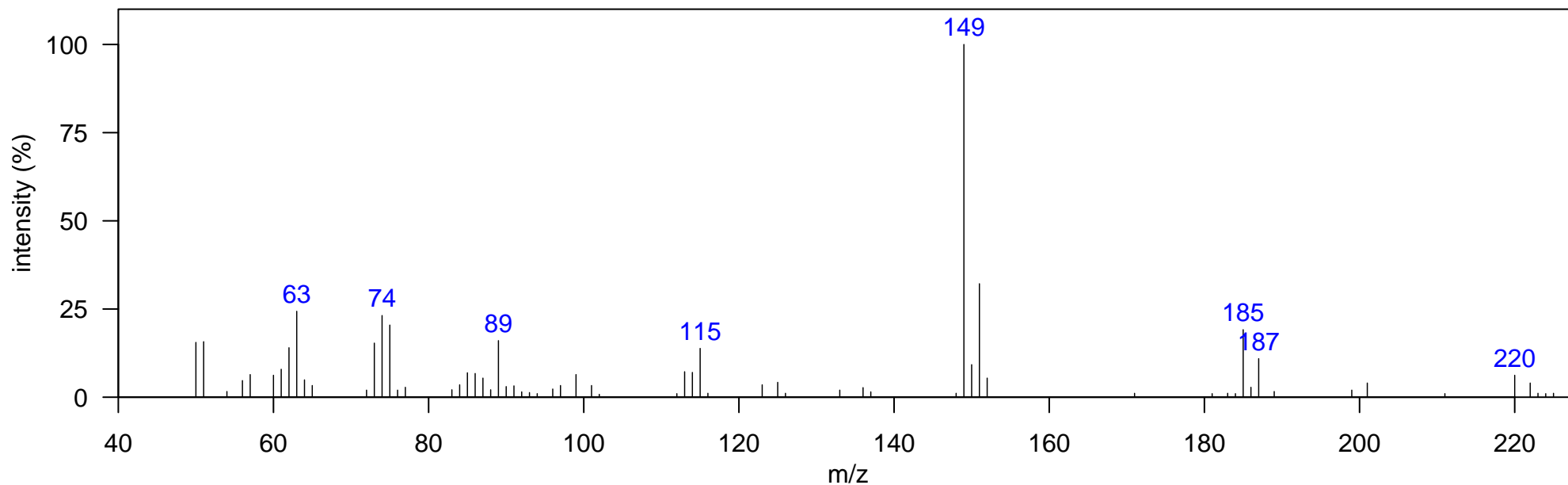

Isomer Information	Isomer				
	InDolphin	Monitored	RT1D	RT2D	
1	trichloro-(propenyl) benzene	FALSE	FALSE	960.294	1.003

m/z [Fragment]
150 [M-Cl <sub>2</sub> ] <sup>+</sup>
185 [M-Cl] <sup>+</sup>
220 M <sup>+</sup>

Name: unknown 1


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

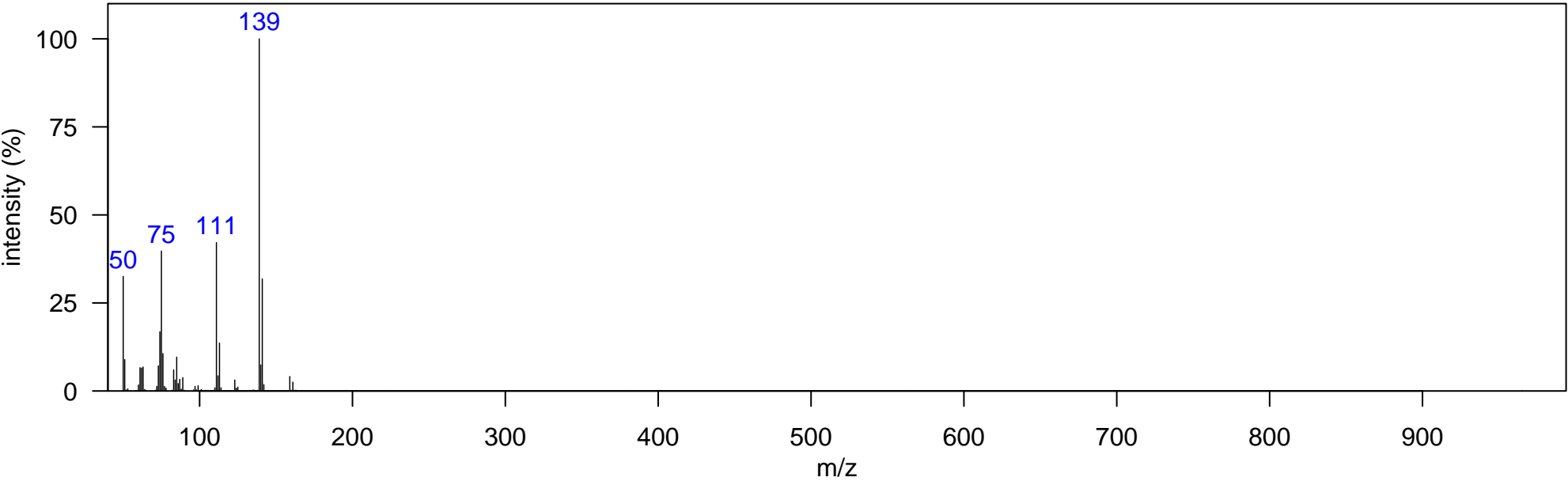

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 1-1	FALSE	FALSE	960.294	1.069
2	unknown 1-2	FALSE	FALSE	970.788	1.056
3	unknown 1-3	FALSE	FALSE	977.784	1.056

m/z [Fragment]

Name: unknown 2


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

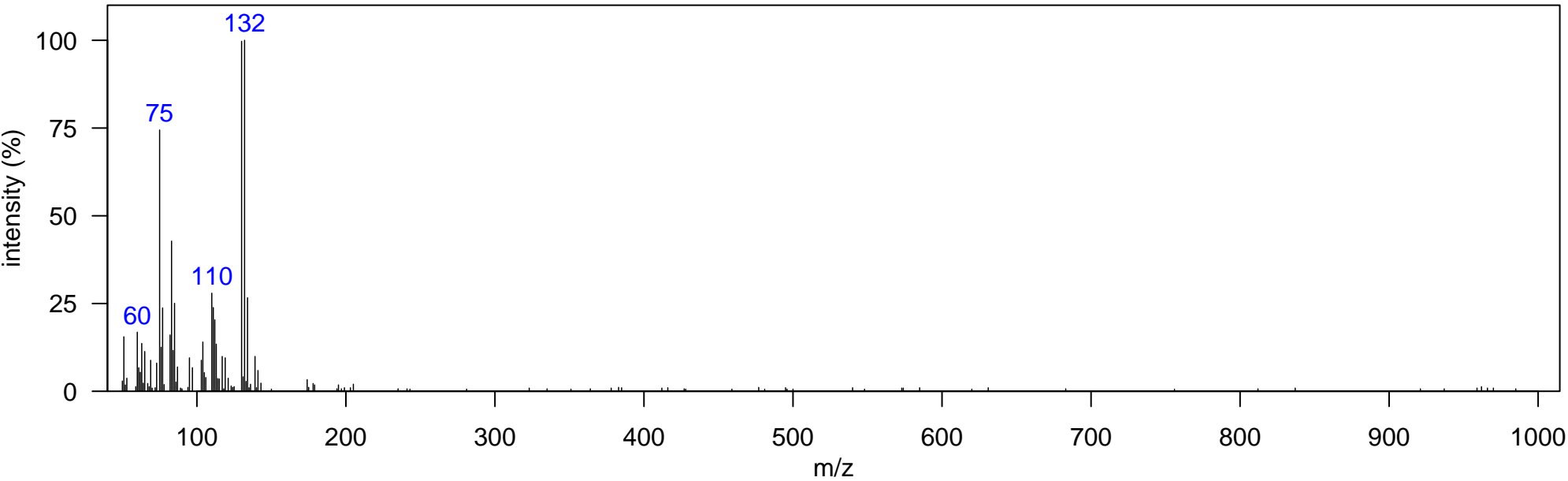

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1 unknown 2	FALSE	FALSE	1026.76	1.01

m/z [Fragment]

Name: pentachloroethylbenzene


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C<sub>8</sub>H<sub>5</sub>Cl<sub>5</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. 1-(Pentachloroethyl)benzene

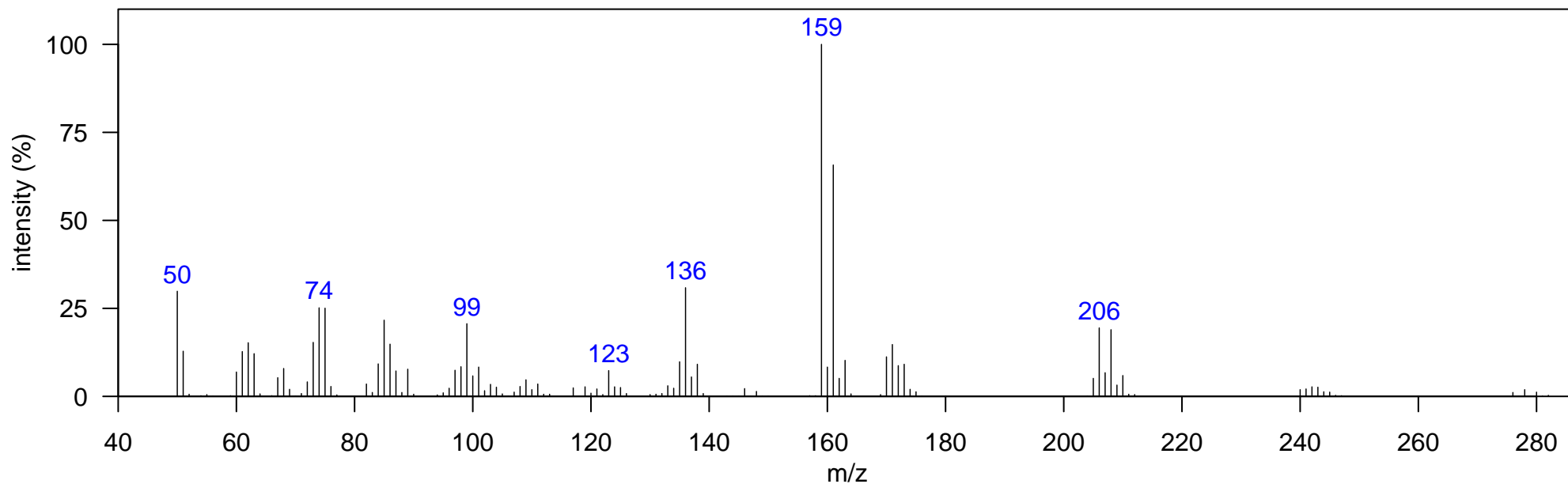

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	pentachloroethylbenzene 1	FALSE	FALSE	1054.74	1.063
2	pentachloroethylbenzene 2	FALSE	FALSE	1093.22	1.082

m/z [Fragment]

159 [M-CCl<sub>4</sub>]<sup>+</sup>  
171 [M-Cl<sub>3</sub>]<sup>+</sup>  
206 [M-Cl<sub>2</sub>]<sup>+</sup>  
241 [M-Cl]<sup>+</sup>  
276 M<sup>+</sup>



Name: unknown 3


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

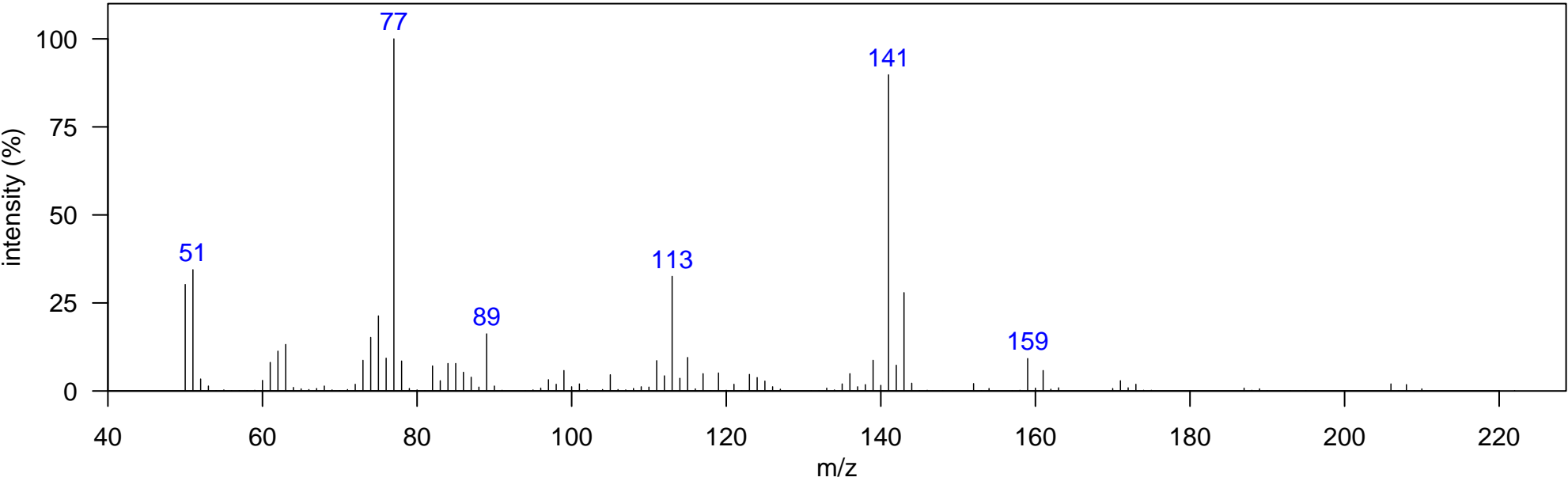

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 3-1	FALSE	FALSE	1086.22	1.129
2	unknown 3-2	FALSE	FALSE	1100.21	1.129
3	unknown 3-3	FALSE	FALSE	1110.71	1.122

m/z [Fragment]

Name: DDM


Technical Mixture: DDT 500PPM

Compound Class: known degradation product

Elemental Formula: C<sub>13</sub>H<sub>10</sub>Cl<sub>2</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. bis(p-chlorophenyl)-methane

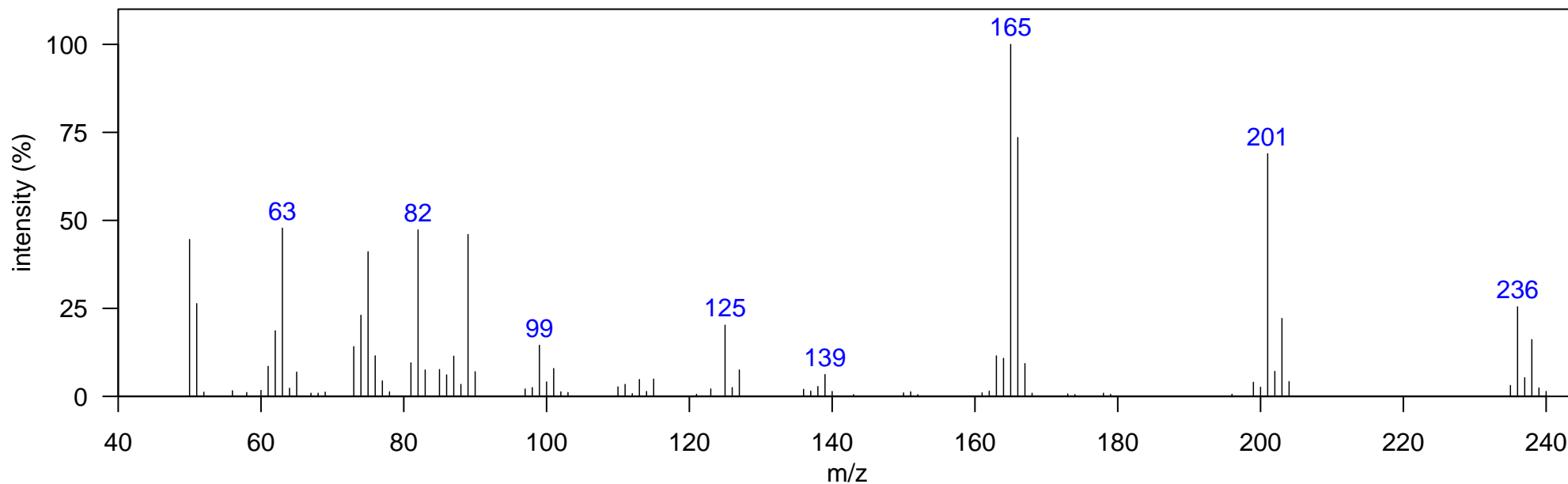

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	DDM 1	TRUE	FALSE	1187.66	1.115
2	p,p'-DDM	TRUE	FALSE	1208.65	1.109

m/z [Fragment]

236 [M]<sup>+</sup>  
201 [M-Cl]<sup>+</sup>  
165 [M-HCl<sub>2</sub>]<sup>+</sup>

Name: unknown 4


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

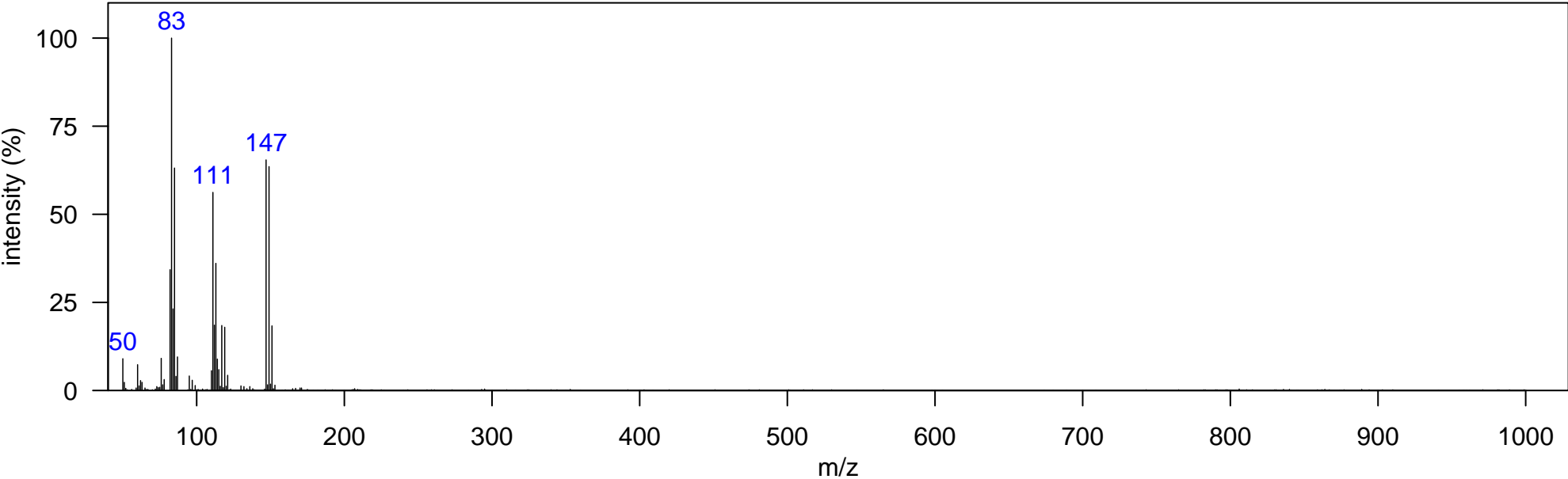

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 4-1	FALSE	FALSE	1219.15	1.102
2	unknown 4-2	FALSE	FALSE	1229.64	1.049
3	unknown 4-3	FALSE	FALSE	1254.13	1.096
4	unknown 4-4	FALSE	FALSE	1551.46	1.320
5	unknown 4-5	FALSE	FALSE	1603.93	1.452

m/z [Fragment]

Name: DDNS


Technical Mixture: DDT 500PPM

Compound Class: known degradation product

Elemental Formula: C<sub>14</sub>H<sub>12</sub>Cl<sub>2</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment: NIST I.D. 1,1-bis(p-chlorophenyl)-ethane; also called DDEt.

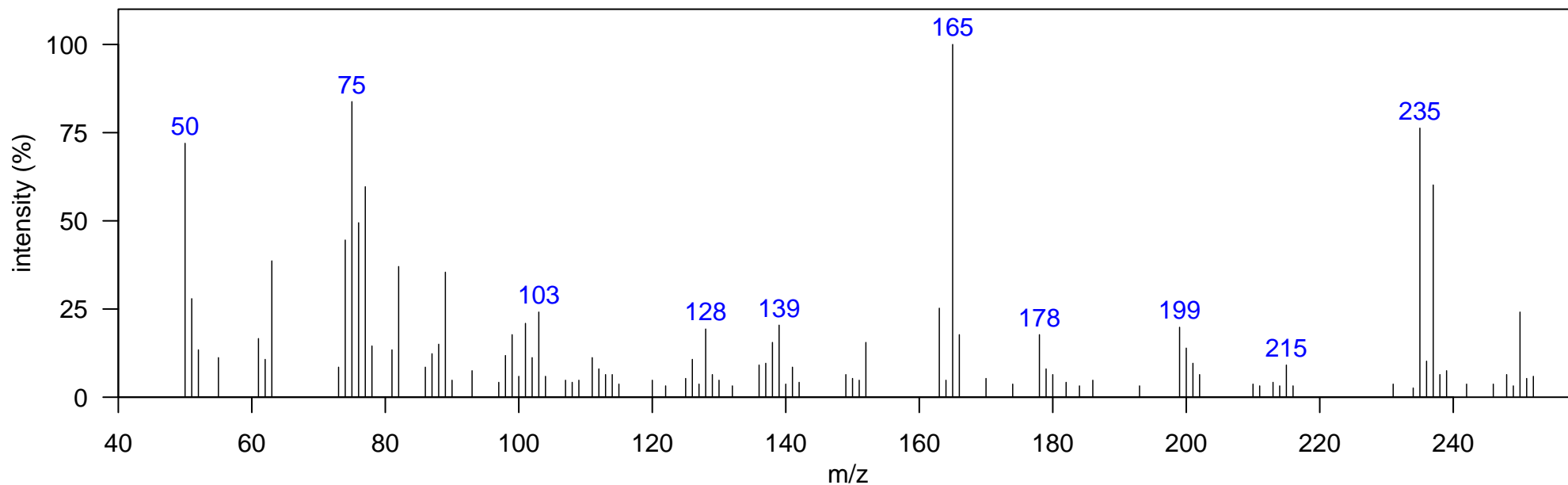

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1	DDNS	TRUE	FALSE	1240.13 1.089

m/z [Fragment]

165 [M-CH<sub>3</sub>Cl<sub>2</sub>]<sup>+</sup>  
235 [M-CH<sub>3</sub>]<sup>+</sup>  
250 M<sup>+</sup>

Name: DDNU


Technical Mixture: DDT 500PPM

Compound Class: known degradation product

Elemental Formula: C<sub>14</sub>H<sub>10</sub>Cl<sub>2</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

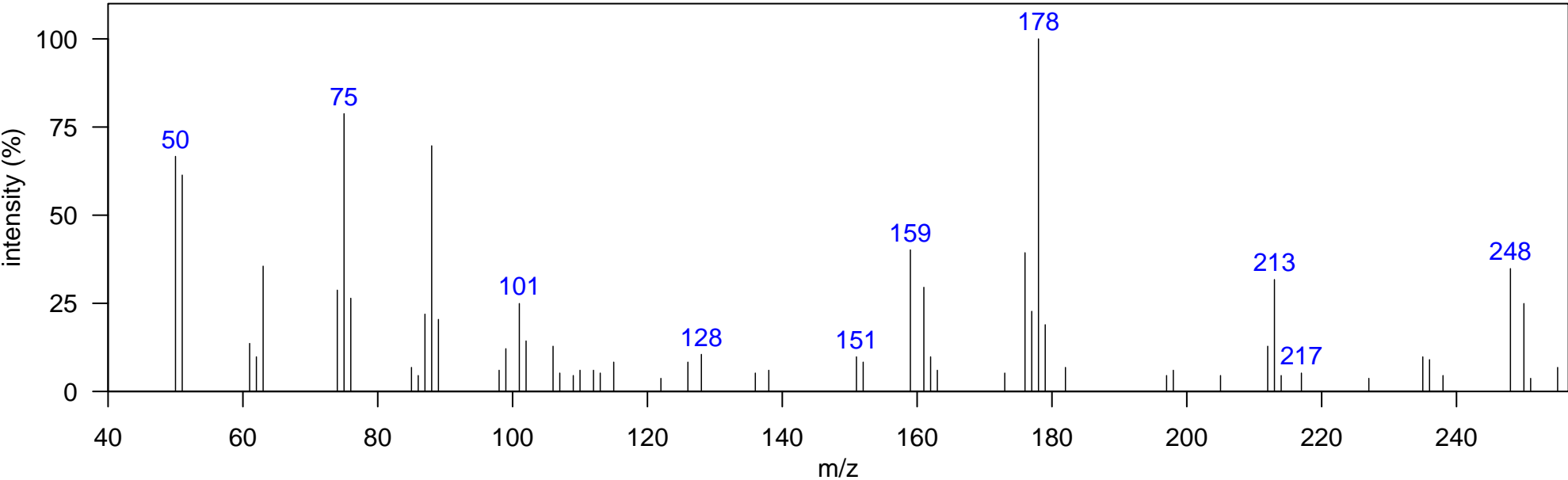

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1 p,p'-DDNU	TRUE	FALSE	1243.63	1.109

m/z [Fragment]

178 [M-Cl<sub>2</sub>]<sup>+</sup>  
213 [M-Cl]<sup>+</sup>  
248 M<sup>+</sup>

Name: DBP


Technical Mixture: DDT 500PPM

Compound Class: known degradation product

Elemental Formula: C13H8Cl2O

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

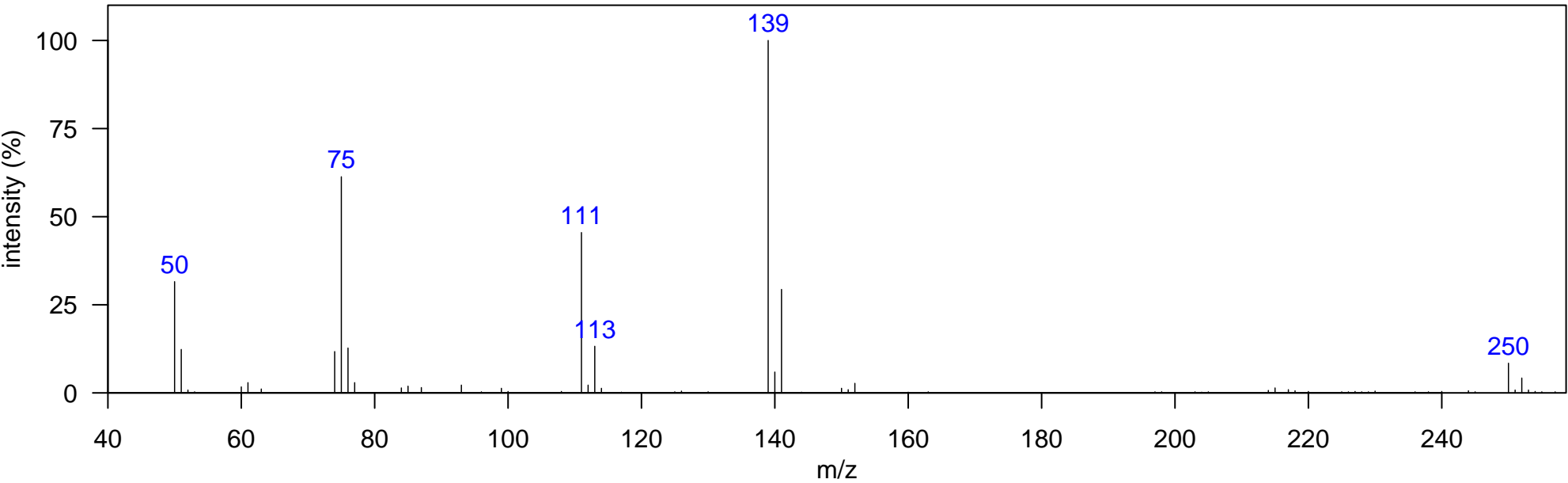


Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	o,p'-DBP	FALSE	FALSE	1254.13	1.181
2	p,p'-DBP	TRUE	FALSE	1285.61	1.162

m/z [Fragment]

75 [M-C7H5Cl2O]+  
111 [M-C7H4ClO]+  
139 [M-C6H4Cl]+  
250 M+

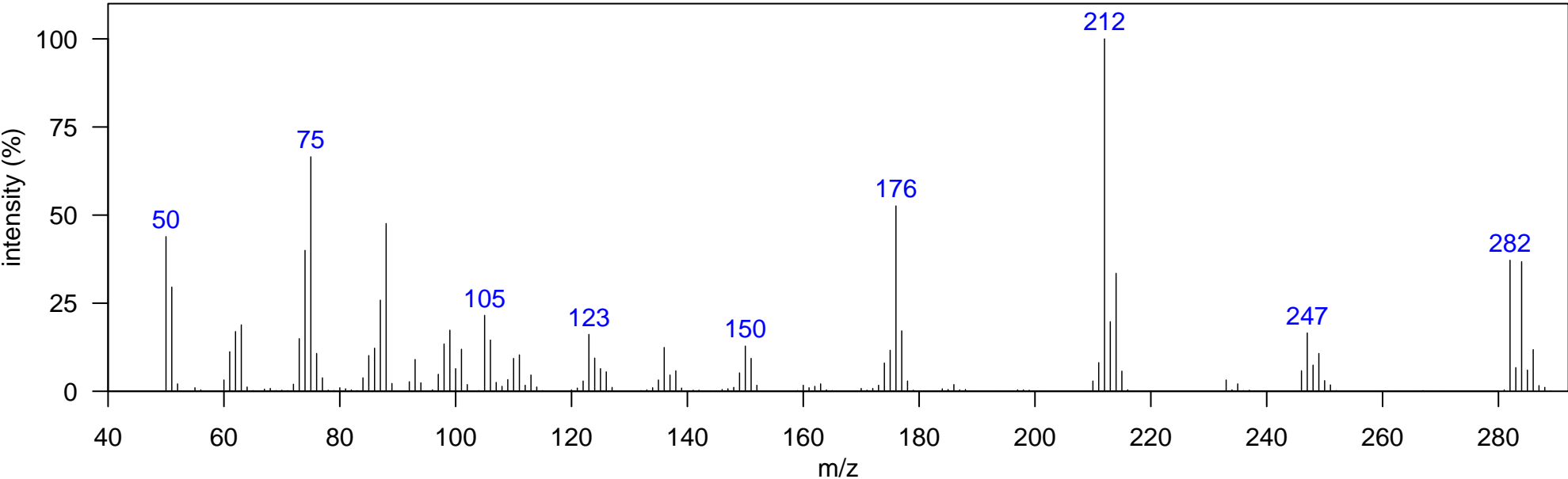

Isomer Information	Isomer InDolphin Monitored				RT1D	RT2D
	1	DDMU 1	FALSE	FALSE	1268.12	1.168
	2	DDMU 2	FALSE	FALSE	1285.61	1.135
	3	DDMU 3	FALSE	FALSE	1296.10	1.162
	4	DDMU 4	FALSE	FALSE	1303.10	1.135
	5	DDMU 5	TRUE	FALSE	1320.59	1.142
	6	DDMU 6	FALSE	FALSE	1331.08	1.195
	7	p,p'-DDMU	TRUE	TRUE	1345.07	1.155
	8	DDMU 7	FALSE	FALSE	1369.56	1.195
	9	DDMU 8	FALSE	FALSE	1376.56	1.195
	10	DDMU 9	FALSE	FALSE	1418.53	1.016
	11	DDMU 10	FALSE	FALSE	1425.53	1.228
	12	DDMU 11	FALSE	FALSE	1436.02	1.016
	13	DDMU 12	FALSE	FALSE	1450.01	1.175
	14	DDMU 13	FALSE	FALSE	1453.51	0.997
	15	DDMU 14	TRUE	FALSE	1467.50	1.155
	16	DDMU 15	FALSE	FALSE	1474.50	1.010

m/z [Fragment]
176 [M-HCl3]+
212 [M-Cl2]+
247 [M-Cl]+
282 M+

Name: unknown 5


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

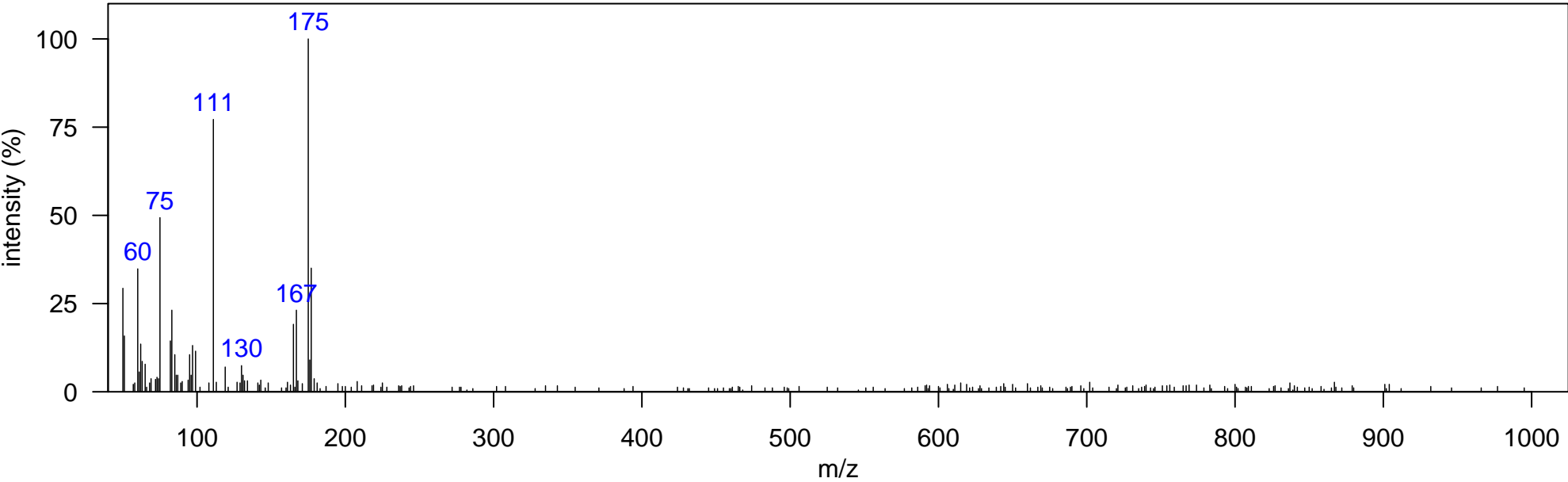

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 5-1	FALSE	FALSE	1292.60	1.162
2	unknown 5-2	FALSE	FALSE	1551.46	1.360

m/z [Fragment]



Name: DDE


Technical Mixture: DDT 500PPM

Compound Class: major degradation product

Elemental Formula: C<sub>14</sub>H<sub>8</sub>Cl<sub>4</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

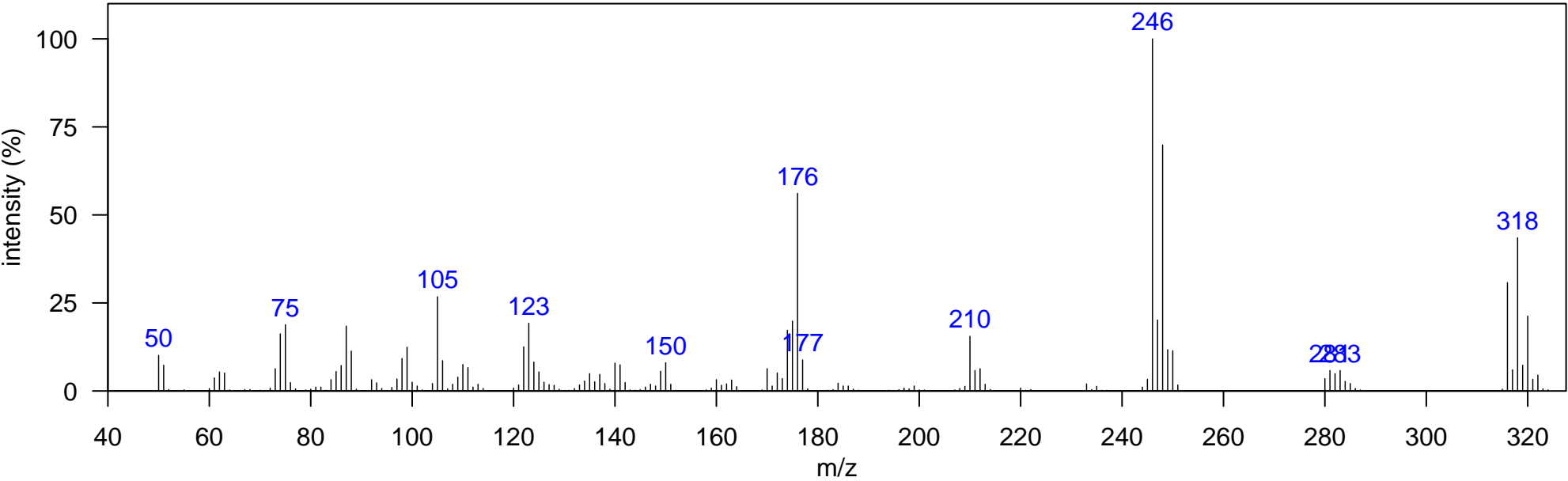

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	DDE 1	FALSE	FALSE	1338.08	1.175
2	o,p'-DDE	TRUE	TRUE	1352.07	1.175
3	DDE 2	FALSE	FALSE	1369.56	1.155
4	p,p'-DDE	TRUE	TRUE	1383.55	1.181
5	DDE 3	FALSE	FALSE	1397.54	1.148

m/z [Fragment]

- 176 [M-Cl<sub>4</sub>]<sup>+</sup>
- 211 [M-Cl<sub>3</sub>]<sup>+</sup>
- 246 [M-Cl<sub>2</sub>]<sup>+</sup>
- 281 [M-Cl]<sup>+</sup>
- 316 M<sup>+</sup>

Name: DDOH


Technical Mixture: DDT 500PPM

Compound Class: known degradation product

Elemental Formula: C14H12Cl2O

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

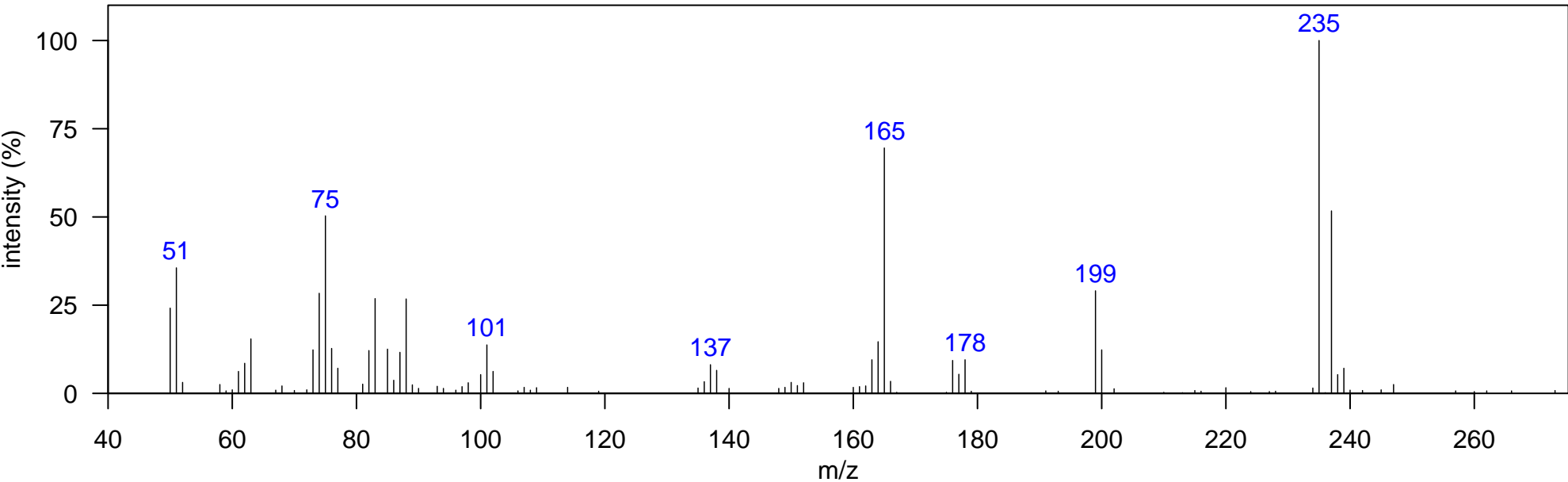

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1	DDOH	FALSE	FALSE	1366.06 1.241

m/z [Fragment]

- 165 [M-CH4Cl2O] +
- 199 [M-CH4ClO] +
- 235 [M-CH3O] +
- 266 M +

Name: unknown 6


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

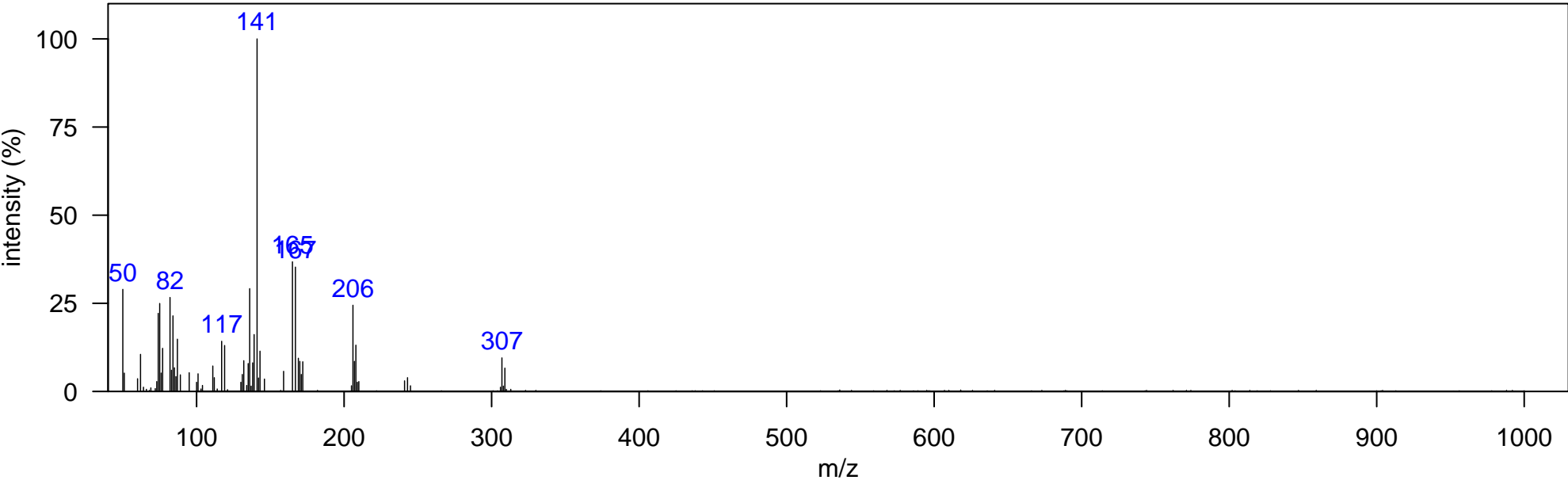

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 6-1	FALSE	FALSE	1380.05	1.162
2	unknown 6-2	FALSE	FALSE	1415.03	1.162

m/z [Fragment]

Name: unknown 7


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

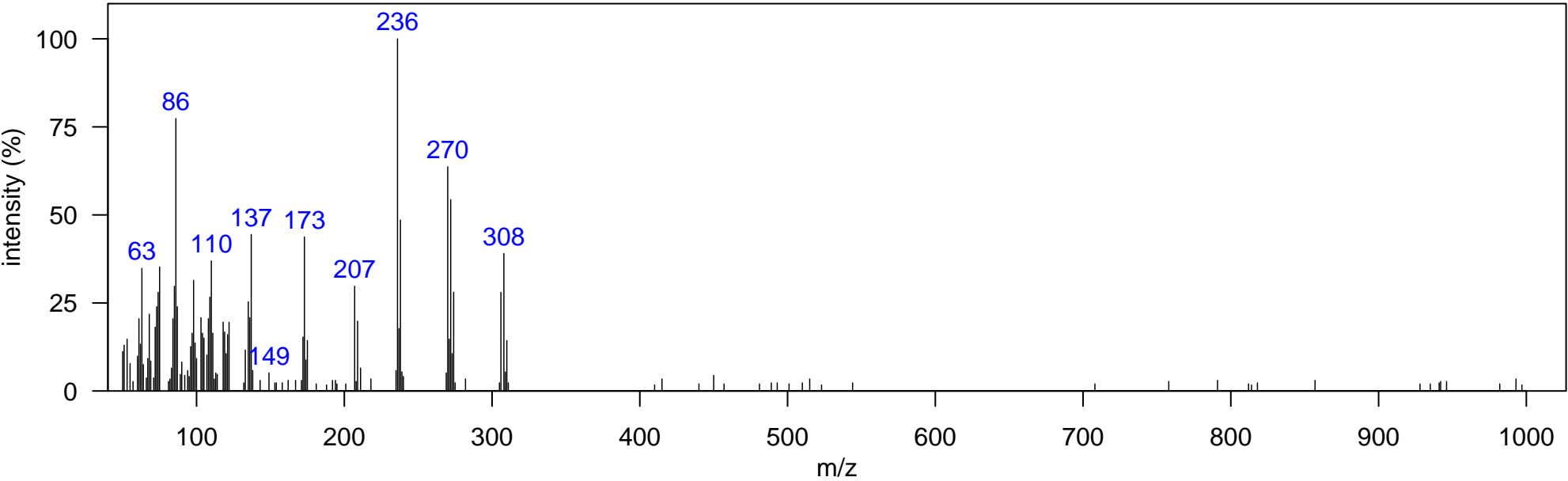

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1 unknown 7	FALSE	FALSE	1390.55	1.261

m/z [Fragment]

Name: DDD


Technical Mixture: DDT 500PPM

Compound Class: major degradation product

Elemental Formula: C<sub>14</sub>H<sub>10</sub>Cl<sub>4</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

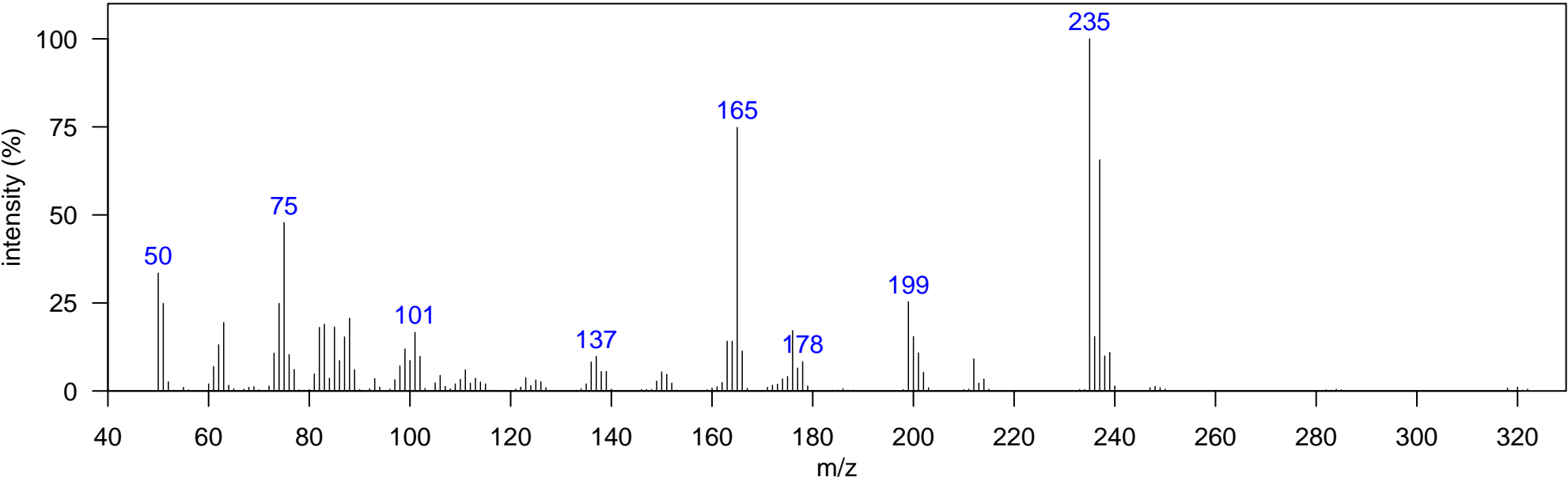

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	o,p'-DDD	TRUE	TRUE	1394.05	1.221
2	DDD 1	FALSE	FALSE	1411.54	1.214
3	p,p'-DDD	TRUE	TRUE	1436.02	1.313
4	DDD 2	FALSE	FALSE	1467.50	1.168

m/z [Fragment]

165 [M-CHCl<sub>4</sub>]<sup>+</sup>  
235 [M-CHCl<sub>2</sub>]<sup>+</sup>  
318 M<sup>+</sup>

Name: unknown 8


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

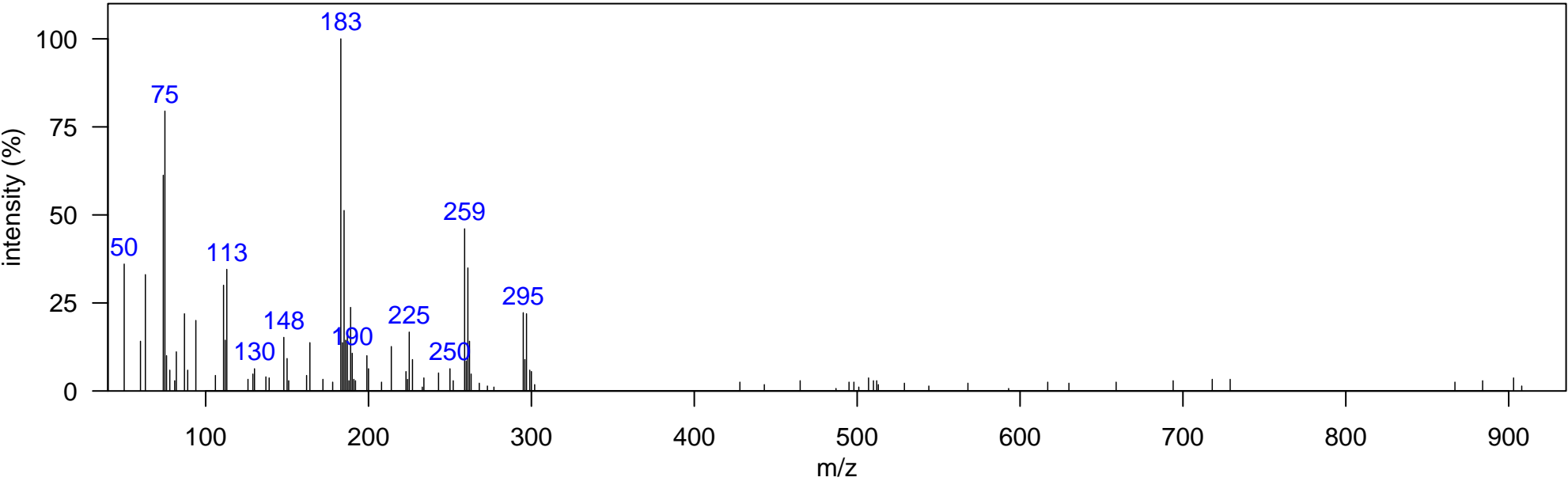

Isomer Information

Isomer	InDolphin	Monitored	RT1D	RT2D
1 unknown 8	FALSE	FALSE	1408.04	1.168

m/z [Fragment]

Name: dichlorodiphenyl


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula: C12H8Cl2O2S

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

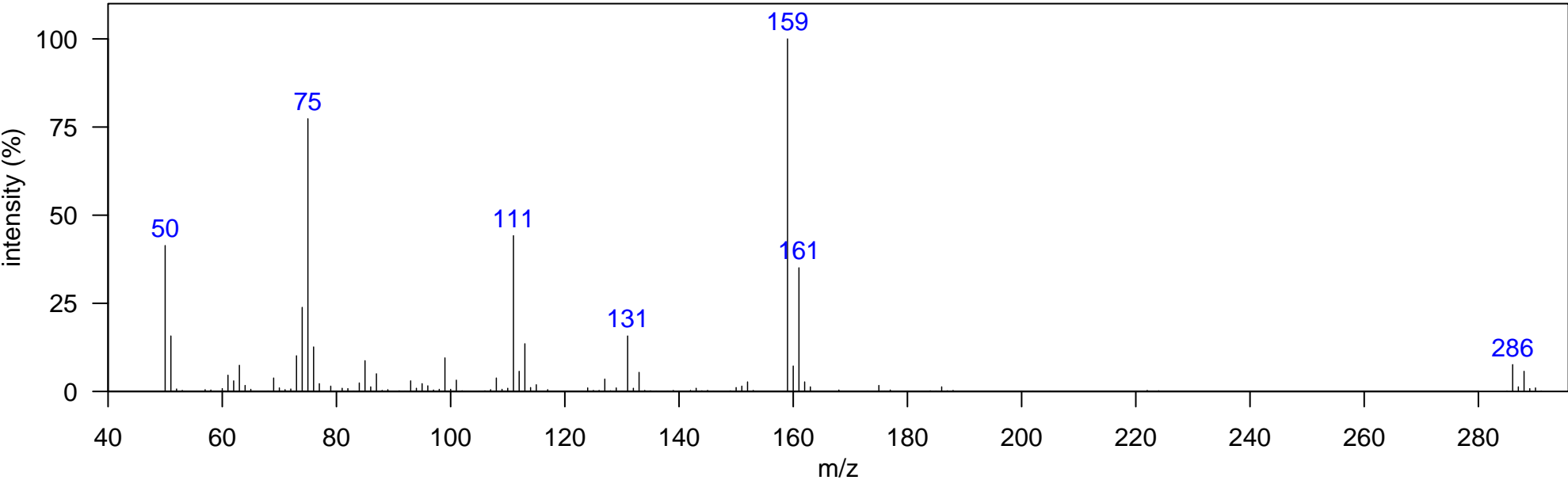

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	p,p'-dichlorodiphenyl sulfone	FALSE	FALSE	1415.03	1.313
2	dichlorodiphenyl sulfone	FALSE	FALSE	1429.03	1.393

m/z [Fragment]

159 [M-C6H4ClSO]+  
286 M+

Name: DDT


Technical Mixture: DDT 500PPM

Compound Class: DDT

Elemental Formula: C<sub>14</sub>H<sub>9</sub>Cl<sub>5</sub>

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

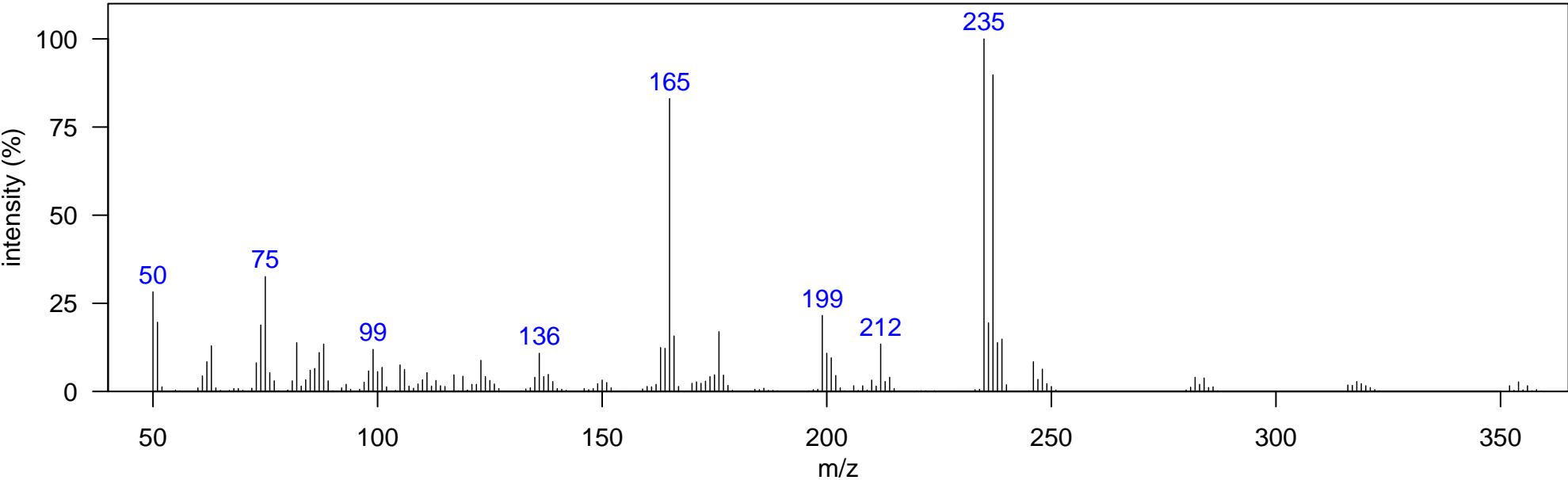

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	DDT 1	FALSE	FALSE	1418.53	1.234
2	o,p'-DDT	TRUE	TRUE	1436.02	1.313
3	DDT 2	FALSE	FALSE	1453.51	1.241
4	p,p'-DDT	TRUE	TRUE	1471.00	1.340

m/z [Fragment]

235 [M-CH<sub>2</sub>Cl<sub>3</sub>]<sup>+</sup>  
354 M<sup>+</sup>



Name: TCPM


Technical Mixture: DDT 500PPM

Compound Class: TCPM

Elemental Formula: C19H13Cl3

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

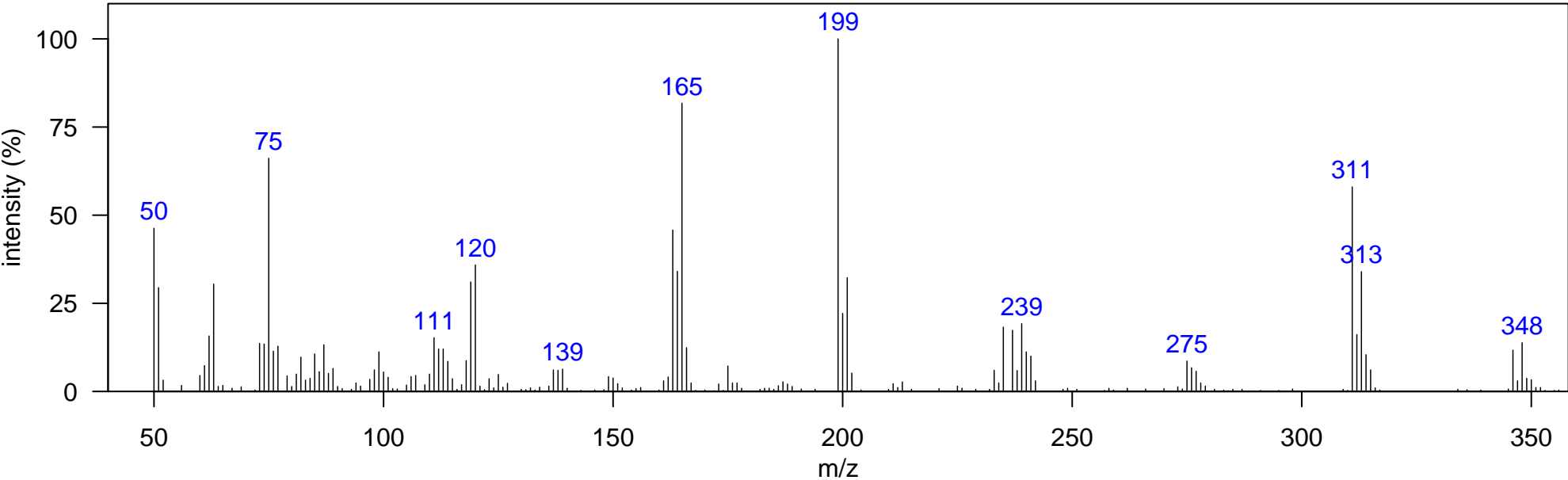

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	TCPM 1	FALSE	FALSE	1540.96	1.426
2	TCPM 2	FALSE	FALSE	1579.44	1.472

m/z [Fragment]

275 [M-HCl2]+  
311 [M-Cl]+  
346 M+

Name: unknown 9


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

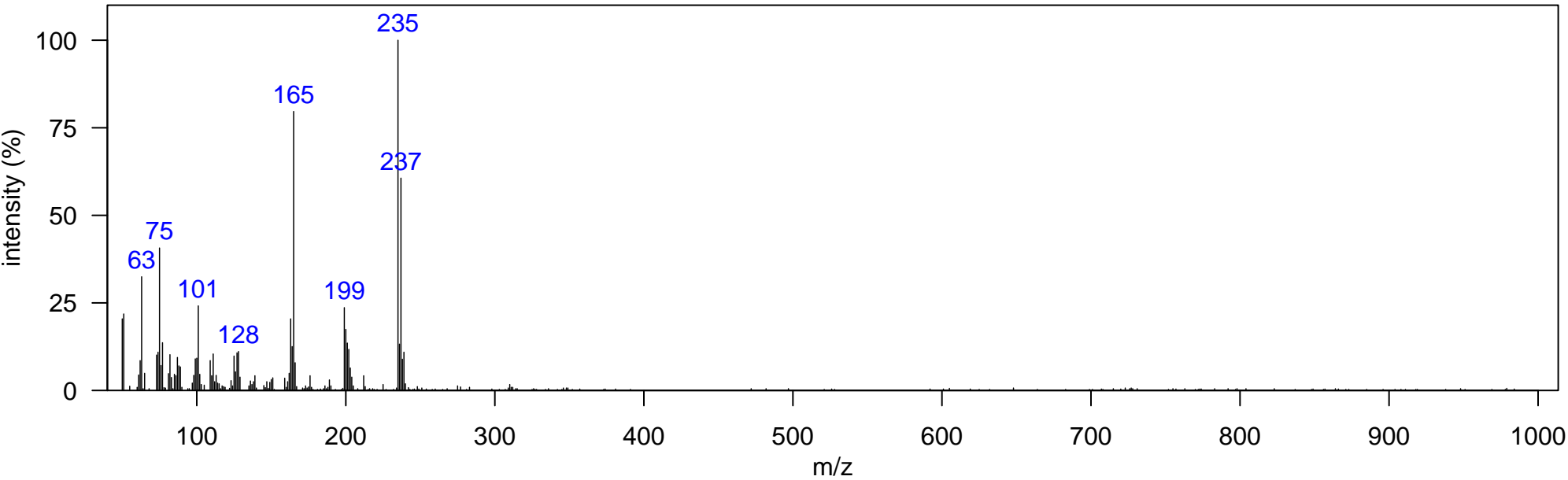

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 9-1	FALSE	FALSE	1537.46	1.360
2	unknown 9-2	FALSE	FALSE	1547.96	1.393
3	unknown 9-3	FALSE	FALSE	1568.95	1.406

m/z [Fragment]

Name: unknown 10


Technical Mixture: DDT 500PPM

Compound Class: other

Elemental Formula:

Instrument: GCxGC-TOF, EI, 70 eV

Comment:

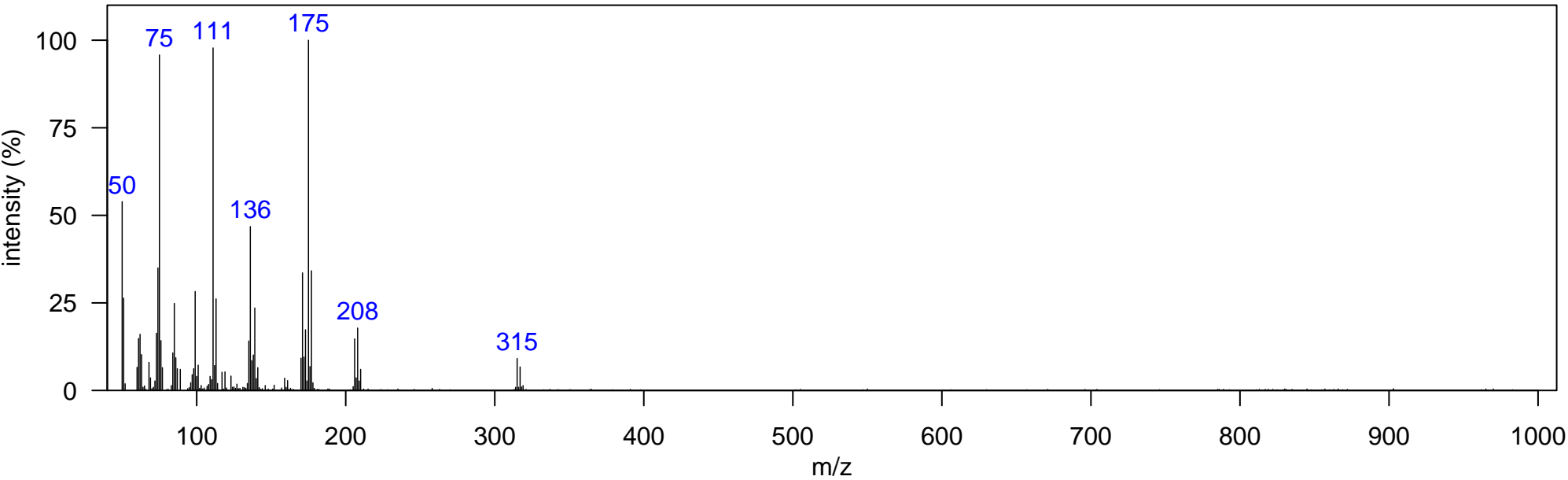

Isomer Information

	Isomer	InDolphin	Monitored	RT1D	RT2D
1	unknown 10-1	FALSE	FALSE	1600.43	1.683
2	unknown 10-2	FALSE	FALSE	1614.42	1.657
3	unknown 10-3	FALSE	FALSE	1635.41	1.927

m/z [Fragment]